MANUAL

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MODE D'EMPLOI

INSTRUCCIONES

Love Magic®

Rechargeable Personal Massager

Model HV-270

HOW TO USE


Follow these simple steps and you will enjoy a Love Magic HV-270 in time.

1. Press the button "⏻" to turn on/off.
2. Press the button "🌊" to enjoy 20 wave pulsations.
3. Press the button "📶" to adjust 8 speeds of continuous vibrations

CHARGING YOUR LOVE MAGIC

1. When your wand's vibration is weak or does not turn on, but the LED is still on. Please recharge.
2. Insert the DC plug into the DC port of the wand.
3. Insert the USB cable into a USB power outlet. The LED on the interface of your Love Magic will light to indicate it is charging.

POSITIONS OF MASSAGING

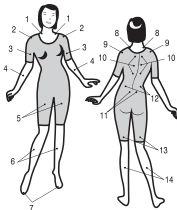Shoulders	1,2,8,9
Arms	3,4
Back muscles	10
Loins	11,12
Legs	5,6,7,13,14

CAUTION: THE MASSAGER SHOULD NOT BE USED OVER SWOLLEN OR INFLAMED AREAS OR SKIN ERUPTIONS. DO NOT USE IN UNEXPLAINED CALF PAIN. CONSULT PHYSICIAN.

IMPORTANT: Any individual who may be pregnant, has a pacemaker, suffers from diabetes, phlebitis, and/or thrombosis, is at an increased risk of developing blood clots due to recent surgery, or who has pins/screws or other medical devices in their body should consult with a physician before using a massaging device designed for home use.

NAME OF PARTS

SPECIFICATIONS

Power supply	Main body	Rechargeable Battery (Lithium-Ion battery 3.7V Typical capacity 2000mAh)
	Charging adapter	Input: AC100-250V 50-60Hz Output: 5V 500mA
Battery charging time	Approx. 3 hours for full charge	
Battery charge duration	Approx. 3 hours	
Weight	Main body: Approx. 460g.	
Dimension	Approx. 13" x 2.4"	
Material	Main body: ABS plastic Head: Silicon	
Charging adapter cord length	Approx. 6 feet (1.8)	
Head vibration level	Approx. 2000r/min,-9000r/min.	
Operating temperature	Approx. 0°C~40°C (32~104°F)	

CARE AFTER USE AND STORAGE, CLEANING

1. Don't store your massager in a moist area, nor in bright sunlight.
2. If necessary, wipe it with a damp cloth certainly not while it's still plugged into the power source. Don't use substances such as benzine, lacquer, thinner or other volatile chemicals. Don't use them on anything electrical.

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following. Read all instructions before using the massager.

DANGER TO REDUCE THE RISK OF ELECTRIC SHOCK:

1. Be extremely carefull never to express the product to any liquid, including body fluids, water, lubricants or cleaners. Serious injury may occur if fluid enter into the product.
The product is NOT waterproof.
2. For external use only. Do not insert into any body cavity.
3. Always unplug this product from the outlet immediately after using and before cleaning.
4. Do not use while bathing or in the shower.
5. Do not place ore store massager where it can fall or be pulled into a tub or sink.
6. Do not place in or drop into water or other liquid.
7. Do not reach for massager that has fallen into water. Unplug immediately.
8. Do not use on any wet body parts. Do not use in the rain.
9. Do not wash the massager with water.
10. Do not modify or dismantle lithium-ion battery.

⚠ WARNING

TO REDUCE THE RISK OF BURNS, FIRE, ELECTRIC SHOCK OR INJURY TO PERSONS:

1. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock or injury to person.
2. Closed supervision is necessary when the massager is used by or near children. This product is for adult use only.
3. Use the massager only for its intended use as described in this manual.
4. Persons who are unable to operate this product by themselves should not use it.
5. Persons who have been told by a doctor to avoid massage should not use this product.
6. Do not use attachments which are not recommended by the manufacturer.
7. Never use this product near a source of heat, such as a heater.
Never use the products when smoking.
8. Keep the cord away from heated surfaces.
9. Never drop or insert any object into any opening.

PRECAUTIONS

1. You'll want to use your massager on your shoulders, arms, back muscles and legs. It's not for your chest and certainly not for use around your thyroid gland (just below the Adam's apple).
2. Don't grasp the vibrating head while your massager is turned ON.
3. The rated maximum continuous use of your massager is 25 minutes.
4. Don't turn the vibrating head by hand or press it tightly to your body. You could bend the head-supporter, and heavy pressure does not produce a stronger effect anyway.
5. Don't use your massager over swollen or inflamed areas, or skin eruptions. Don't use on an unexplained calf pain. Whenever in doubt, consult your physician.
6. After use, store the product in a place out of reach of children.

THE MAINTENANCE FOR LITHIUM-ION BATTERY

1. Recharge the battery as soon as it runs out.
2. Make sure the battery would be charged for at least 1 hour for every three months, because the battery could be dead after a long period self discharge.

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